

Cementation Solutions

RelyX™ Luting 2

Resin Modified Glass Ionomer Cement



Easy just
got easier



3M ESPE

RelyX™ Luting 2 – comfort plus performance

Now even easier to use with
NEW tack light-cure clean-up feature.

RelyX™ Luting 2 Resin Modified Glass Ionomer (RMGI) Cement from 3M™ ESPE™ – well suited for the cementation of PFM restorations and of restorations to implant abutments – offers the perfect combination of comfortable handling and proven clinical performance.

Its advanced fluoride-releasing formula includes initiators that now allow for easy clean-up following a short five second tack light-cure.

Benefit from the convenient Clicker™ handling.

RelyX Luting 2 cement in the handy Clicker™ Dispenser was designed to make your procedure faster and easier. The two different coloured pastes allow for easier monitoring of a homogeneous mix.




Convenient delivery from the Clicker™ Dispenser for 1:1 dosing and easy mixing.

Indications:

RelyX™ Luting 2 Cement is recommended for the cementation of/to:

- PFM crowns and bridges
- Metal crowns, bridges, inlays and onlays
- Endodontic posts
- Orthodontic appliances
- All-zirconia or all-alumina strengthened core-ceramic restorations (such as Lava™, Procera®, AllCeram)
- **NEW!** Implant abutments

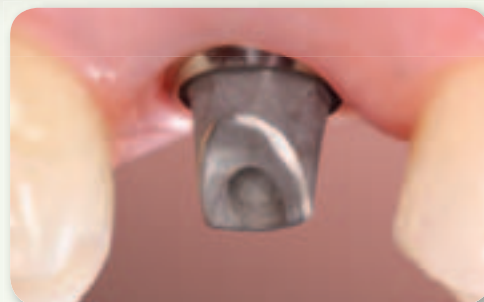

RelyX™ Luting 2 Cement is recommended for the cementation of restorations to implant abutments.



The first RMGI cement with tack light-cure option

Enjoy easy clean-up of excess cement within seconds.


You can speed up the gel state by light-curing excess RelyX™ Luting 2 Cement for about 5 seconds per surface with the Elipar™ S10 Curing Light from 3M ESPE or any other conventional light-curing device. Alternatively, you can wait for the excess to self-cure – about two minutes after seating in the mouth. Once excess RelyX Luting 2 cement reaches a gel state, it can be very easily dislodged in large sections from the margins.




Tack curing of RelyX™ Luting 2 Cement to remove excess easily.



Easy removal of excess cement in one piece.

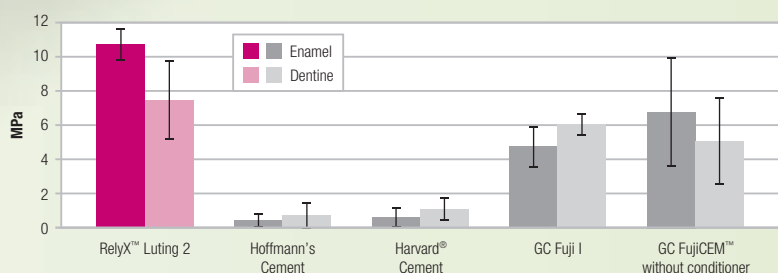

Compared to competitive products, RelyX™ Luting 2 Cement shows high shear bond strength values to enamel/dentine.

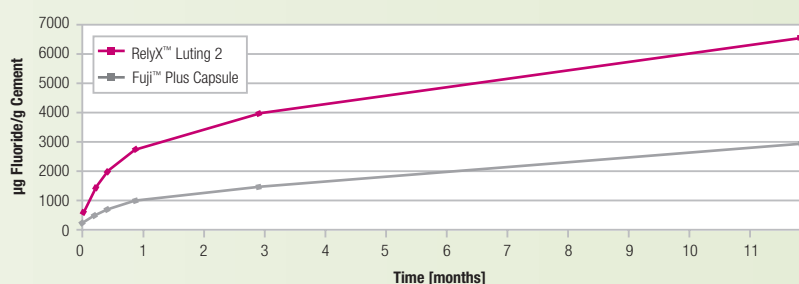

RelyX™ Luting 2 Cement offers sustained fluoride release.

Proven to have better bond strengths.

RelyX Luting 2 cement provides molecular bonding to the tooth surface. Compared to other cements, it delivers exceptional shear bond strength values – without a separate tooth conditioner. And there is virtually no post-operative sensitivity.

Source: 3M ESPE internal data 2012 (Notched-edge adhesion).

Fluoride release.

For your caries-prone patients, RelyX Luting 2 cement offers sustained fluoride release.

Source: Partially cited from: Ton T, Wang Y, Falsafi A, Oxman J. Fluoride Release of a Novel Automix Resin-Modified Glass-Ionomer Luting Cement. J Dent Res 91 (Spec Iss A) 1319, 2012. (www.dentalresearch.org).

Recommendation by Cement & Indication

RelyX™ Luting 2 Clicker™
Resin Modified Glass Ionomer Cement

RelyX™ Ultimate Clicker™
Adhesive Resin Cement

RelyX™ U200
Self-Adhesive Resin Cement

Metal/Metal based			
Inlays/Onlays	+	+	++
Crowns/Bridges	++	+	++
Endodontic Posts	++	+	++
Maryland Bridges	-	++	+
Glass Ceramics (incl. Li. Disilicate)			
Inlays/Onlays	-	++	+
Partial Crowns/Table tops	-	++	+
Crowns/Bridges	-	++	++
Veneers	-	+	-
Oxide Ceramics (e.g. Lava™ Zirconia)			
Crowns/Bridges	+	+	++
Endodontic Posts	+	+	++
Maryland Bridges	-	++	+
Resin Nano Ceramics (e.g. Lava™ Ultimate CAD/CAM Restorative)			
Inlays/Onlays	-	++	-
Crowns	-	++	-
Veneers	-	+	-
Resin Composites			
Inlays/Onlays	-	++	+
Crowns	-	++	++
Endo Posts (e.g. RelyX™ Fiber Post)	-	+	++
Veneers	-	+	-

++ Highly recommended*

+ Recommended

- Not indicated

*Either better performance for this indication or easier handling at equal performance.

Ordering Information



RelyX™ Luting 2 Trial Kit (3525TKA)



RelyX™ Luting Introductory Kit (3505)

Item No. Product Information

RelyX™ Luting 2

3525TKA Trial Kit

1 Clicker™ Dispenser 11 g, 1 mixing pad

3525A Refill

2 Clicker™ Dispenser 11 g each, 1 mixing pad

RelyX™ Luting

3505 Introductory Kit

1 powder 16 g, 1 liquid 9 ml, 1 measuring spoon, 2 mixing pads

3505 P Refill

1 powder 16 g

3505 L Refill

1 liquid 9 ml

3M ESPE

3M Deutschland GmbH
Location Seefeld
3M ESPE · ESPE Platz
82229 Seefeld · Germany
info3mespe@mmm.com
www.3MESPE.com

www.3MESPE.com

3M, ESPE, Clicker, Elipar, Lava and RelyX are trademarks of 3M Company or 3M Deutschland GmbH. Used under license in Canada.

All other trademarks are owned by other companies.

© 2012, 3M. All rights reserved.

XX (0X.2012)